



CONSOLOR LTD
*improving posture
 encouraging comfort*


www.consolor.co.uk

Consolor Shoulder Harness Mounting

The **Consolor Shoulder Harness Mounting** is intended for use as a **postural support positioning device**. It makes it possible to use a **postural shoulder harness** with a back support with **scapula cut outs** or where the back support is shaped or positioned below the top of the shoulders. It is mounted to the vertical headrest tube and adjusted to the desired position, which would usually be level with the top of the shoulders.

The **Consolor Shoulder Harness Mounting** is designed to be used in conjunction with wheelchair back supports fitted with headrests which use a 1/2" square vertical tube.

Fitted with two cam buckles for fixing 1" wide harness webbing straps, these can be adjusted laterally in order that the top straps of the harness are positioned correctly over the client's shoulders. The cam buckles can be replaced with your own fittings if preferred.

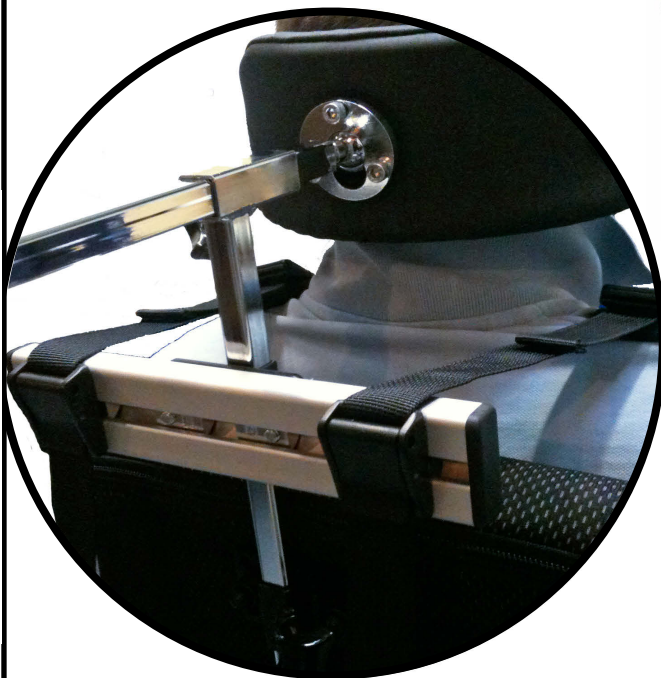

Product codes

CS-SHM-S
 CS-SHM-L

Description

6" wide - ideal for MaTRx backs between 12" and 15" wide
 8" wide - ideal for MaTRx backs between 16" and 20" wide